

Cobolt FIC-01

Mounting plate for laser-to-fiber coupling with beam access

Cobolt model number	FIC-01
Material	Aluminium
FIC-01 is suitable for use with all Cobolt lasers in the following model:	Cobolt 04-01 Series
Dimensions [mm] [inches]	185 x 140 x 51 7.3 x 5.5 x 2.0

Fiber coupling options are available for all Cobolt lasers. The option FIC-01 is recommended for the Cobolt 04-01 Series lasers. The option consists of a mounting plate for the laser and a coupling block to mount the fiber coupler onto. Fiber coupling options for either multimode fiber (MM) or singlemode polarization maintaining (SM/PM) fiber can be ordered separately. The FIC-01 is designed to be bolted to a standard metric or imperial optical table.

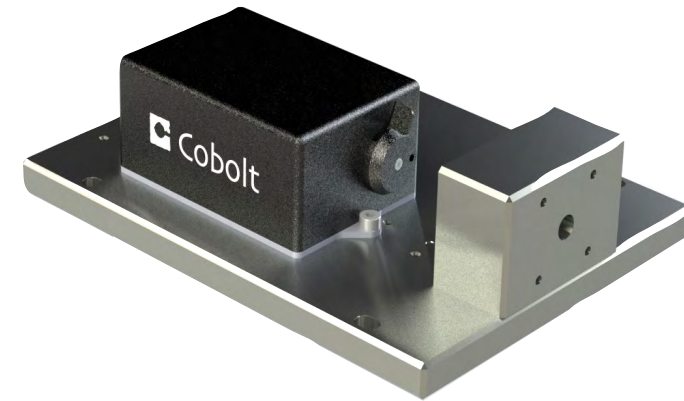


Illustration of possible configuration. Laser and fiber coupling option should be ordered separately.

